
Heraeus

PALACOS® MV+G

Medium viscosity, radiopaque bone cement
containing Gentamicin

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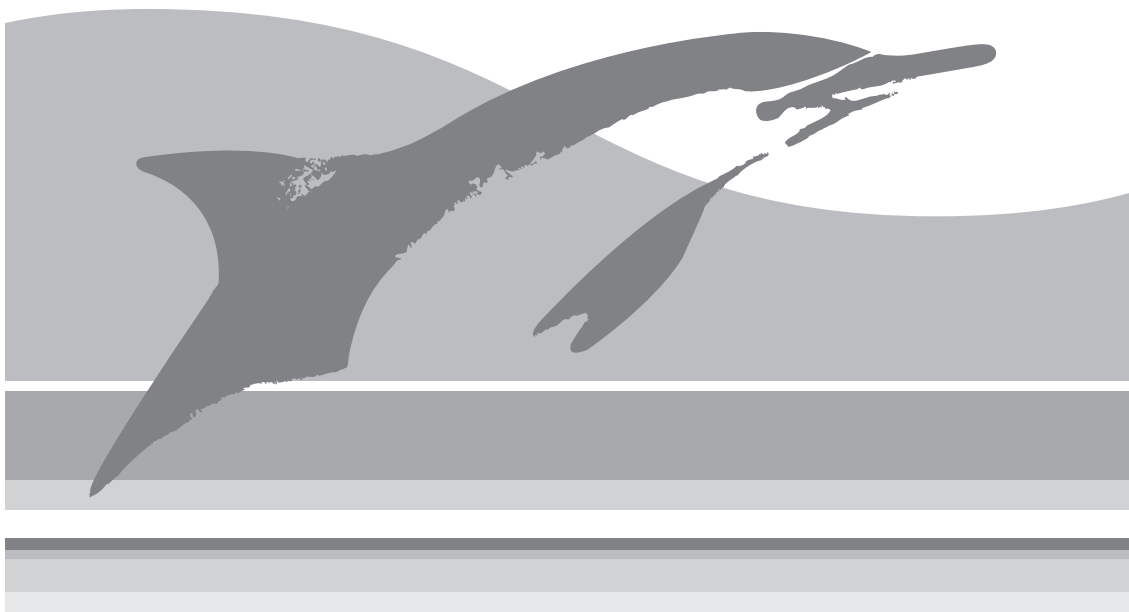PALACOS®MV+G

Properties – Composition

PALACOS®MV+G is a radiopaque, quick-setting bone cement with the addition of **gentamicin** as an antibiotic and exhibits low initial viscosity. It is obtained by mixing a polymer powder component with a liquid monomer component. Zirconium oxide has been added to the cement powder as an X-ray contrast medium. The sterile-filtrated monomer component is supplied in an amber glass ampoule and comes in a sterile blister pack.

The polymer powder component is supplied in a double sterile packaging. The inner polyethylene sachet which contains the powder component is wrapped in an additional polyethylene sachet; both sachets were sterilized with ethylene oxide. The polyethylene sachets are contained in a non-sterile protective aluminium packaging. The polymer powder component is supplied in a double sterile packaging. The inner polyethylene sachet which contains the powder component is wrapped in an additional polyethylene sachet; both sachets were sterilized with ethylene oxide. The polyethylene sachets are contained in a non-sterile protective aluminium packaging. Chlorophyll (Type E141) has been used to obtain the green colour of **PALACOS®MV+G** in order to ensure clear visibility of the cement at the operating site.

After mixing, a plastic dough is obtained which is filled into the bone as an anchoring medium. The cement which then hardens in the bone allows stable fixation of the endoprostheses. The stress forces resulting from motions are transferred via the cement coating widely onto the bone.

Intended use

PALACOS®MV+G is a radiopaque cement-like substance which allows the implantation and fixation of prostheses in the bone.

Indications

PALACOS®MV+G is indicated for the fixation of prostheses in the bone in partial and total arthroplastic surgery of the hip, knee or other joints if an infection with gentamicin-sensitive germs is present or suspected.

PALACOS®MV+G offers protection against accumulation of gentamicin-sensitive germs on the graft and the adjacent tissue.

Contraindications

PALACOS®MV+G must not be used during pregnancy or nursing. In cases of known hypersensitivity to the constituents of the bone cement **PALACOS®MV+G** must not be used.

PALACOS®MV+G must not be used in cases of serious renal insufficiency.

Warning information – Side effects

PALACOS®MV+G has not been evaluated with regard to spinal surgery. In some cases, the use of this cement beyond the listed indications in spine surgery resulted in serious, life-threatening complications. Cases of pulmonary embolism, respiratory and cardiac insufficiency and death have been reported.

Prior to using **PALACOS®MV+G** the surgeon should be familiar with its properties, handling and application during arthroplastic surgery. It is also recommended for surgeons to practice mixing, handling and application of **PALACOS®MV+G** prior to the use. Precise knowledge is also required, if mixing systems and syringes are used for the application of the cement. The monomer liquid is highly volatile and flammable; accordingly, suitable precautionary measures should be taken for use in the operating room. The monomer is also a powerful lipid

Composition	PALACOS®MV+G 1x40
	1 sachet of 44,9 g powder contains:
Poly (methyl acrylate, methyl methacrylate)	38.3 g
Zirconium dioxide	5.3 g
Benzoyl peroxide	0.4 g
Gentamicin base (as sulphate)	0.55 g
Methyl methacrylate	18.4 g
N,N-Dimethyl-p-toluidine	0.4 g
Other constituents:	In the powder: colourant E141 In the liquid: colourant E141, hydroquinone

solvent and should not come into direct contact with the body. When working with the monomer or the cement, gloves must be worn to ensure adequate protection against the penetration of the monomer (methyl methacrylate) into the skin. PVP (three-layer polyethylene, ethylene-vinyl alcohol-copolymer, polyethylene) and Viton/Butyl gloves have proven to provide good protection over an extended period. Wearing two pairs of gloves has also proved to offer adequate protection.

The use of latex or polystyrene-butadiene gloves, however, must be avoided. Please request the confirmation of your glove supplier whether the respective gloves are suitable for use with **PALACOS®MV+G**.

The monomer vapors may irritate the respiratory tract and the eyes and possibly damage the liver. Skin irritations have been reported, which must be attributed to the contact with the monomer.

Manufacturers of soft contact lenses recommend the removal of the lenses in the presence of harmful or irritating vapors. Since contact lenses are permeable to liquids and gases, they should not be worn in the operating room if methyl methacrylate is used.

Precautionary measures

Blood pressure, pulse and breathing must be carefully monitored during and immediately after implanting the bone cement. Any significant change of these vital signs must be immediately responded to with adequate measures.

If **PALACOS®MV+G** is used for a total hip endoprosthesis, the proximal part of the medullary (bone marrow) canal of the femur and the acetabulum need to be thoroughly cleaned, aspirated and dried.

To reduce the considerable increase in pressure in the intraosseous area during the implantation of the prosthesis, it is recommended to use suction drainage. If pulmonary, cardiovascular complications arise, monitoring and – in some cases – even increasing the blood volume may be required. In case of acute respiratory insufficiency, anaesthesiological measures should be taken.

Undesired effects

Frequently, a temporary drop in blood pressure immediately after the implantation of the bone cement and the endoprosthesis has been observed. Rare cases of hypotension, including anaphylactic shock, followed by cardiac arrest and sudden death have been reported.

The following, additional undesired effects of the use of methyl methacrylate bone cement have been observed: Thrombophlebitis, superficial wound infection, deep wound infection, pulmonary embolism, haemorrhage and haematoma, trochanter bursitis, loosening or displacement of the prosthesis, trochanter detachment.

Other side effects observed: heterotopic bone regeneration, myocardial infarction, temporary cardiac arrhythmia, cerebrovascular accident.

When using **gentamicin**, the typical side effect of this antibiotic may occur, in particular damage to hearing and renal damage. The occurrence of this side effect, however, is unlikely due to the very low serum level.

Interactions

When administering muscle relaxants and ether, the neuromuscular blocking properties of the gentamicin may be reinforced; the occurrence of this side effect, however, is unlikely due to the very low serum level.

Incompatibilities

Aqueous solutions (e.g. containing antibiotics) must not be added to the bone cement since they considerably impair the mechanical properties of the cement.

Dosing and preparation

A single dose is prepared by mixing the entire content of a powder sachet with one ampoule. The quantity to be used depends on the respective surgical operation and the technique employed. Before the beginning of the operation, at least one additional dose of **PALACOS®MV+G** should be readily available. Each dose is prepared separately.

The following is required to prepare the bone cement:

Sterile working surface, porcelain or stainless steel bowls, sterile mixing spoons or porcelain or stainless steel spatulas or a sterile vacuum mixing system.

The protective aluminium packaging, the outer non-sterile polyethylene sachet and the blister pack of the ampoules should be opened by an assistant in a way to maintain sterility. The sterile polyethylene sachet and the ampoule are placed under aseptic conditions on a sterile table. Sterile conditions must be ensured when opening the polyethylene sachet and the ampoule.

Application

Two different methods can be used for mixing:

Manual mixing

Mixing under vacuum

Manual mixing

The liquid is poured into a bowl and the powder is added. Do not open the ampoule over the mixing device to prevent contamination of the cement with glass fragments. Then the mixture is stirred carefully for 30 sec.

If the dough-like mass no longer sticks to the rubber gloves, it can be processed. The application time depends on the material and room temperature. If the required consistency is obtained, the cement can be applied. To ensure adequate fixation, the prosthesis should be implanted and stabilized within the given time until the bone cement has hardened completely. Surplus cement must be removed as long as it is still soft.

If additional cement is needed during the operation, another sachet of powder can be mixed with one ampoule of liquid as described above. The kneadable mass which is obtained must be applied to the previously applied cement before it hardens. It is always required to mix the entire content of a sachet with the entire content of an ampoule.

Mixing under vacuum

To obtain a bone cement with reduced porosity, the components are mixed under vacuum. For this purpose an airtight closed system and rapid build-up of sufficient vacuum in the mixing equipment are required (absolute pressure: approx. 200 mbar). The stirring times for mixing under vacuum and mixing without vacuum are identical

(30 sec). For details of the mixing method refer to the instructions of the mixing system used.

Storage

Do not store above 25 °C (77 °F).

Shelf life/Sterility

The shelf life is 3 years. The expiration date is printed on the folded box, the protective aluminium packaging and the inner sachet. **PALACOS®MV+G** must not be used after the expiration date. The contents of unused but opened or damaged packs must be discarded and must not be resterilized. **PALACOS®MV+G** has been sterilized with ethylene oxide gas and must not be resterilized. The polymer powder must not be used if it exhibits yellow discoloration.


Manufacturer: Heraeus Medical GmbH
Registered office address: Philipp-Reis-Straße 8/13,
61273 Wehrheim, Germany
Manufacturing address: Philipp-Reis-Straße 8/13,
61273 Wehrheim, Germany
Telephone: +49 (0) 6181.35-3399
Fax: +49 (0) 6181.35-3366
Website: www.heraeus-medical.com

After sale-service agent:
LandMover Medical & Technology (Beijing) Co., Ltd.
Address: Room 606, 6th Floor, Building 2, No. 106,
Majiapu east road, Fengtai District, Beijing, China
Tel: +86 10 5803 1636
Fax: +86 10 5803 1636 EXT 804


License number: CFDA (I)20113651064 (更)
Product Standard: YZB/GER 3880-2010

Working times for manual mixing (not pre-chilled bone cement)

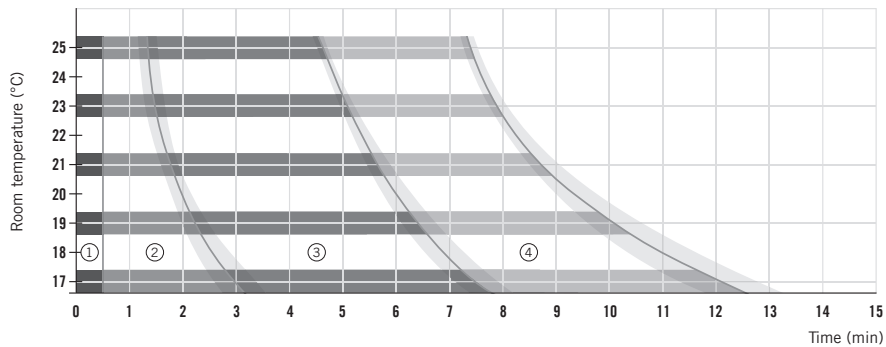Working times for vacuum mixing (not pre-chilled bone cement)

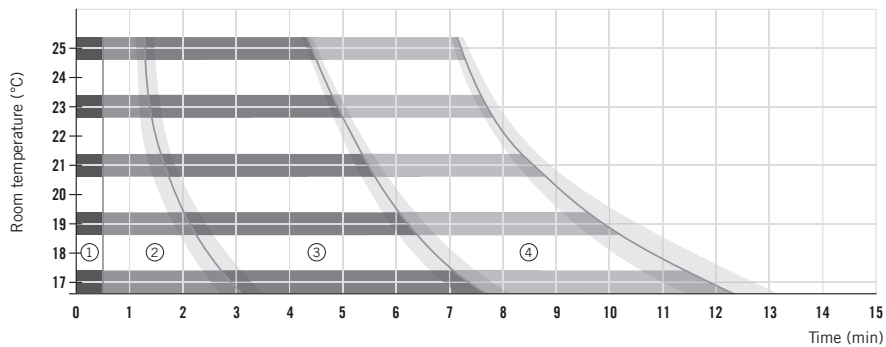■ Mixing ① ■ Waiting ② ■ Application ③ ■ Setting ④

Test conditions: Not pre-chilled vacuum mixing system PALAMIX®, 55% humidity.

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PALACOS® is a trademark of Heraeus Medical GmbH.